



Insert filter holder KL 200

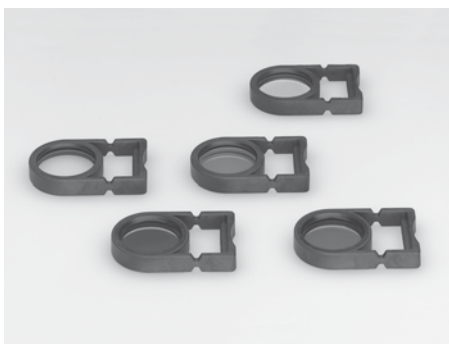

Insert filters for KL 1500 LCD and KL 1500 compact

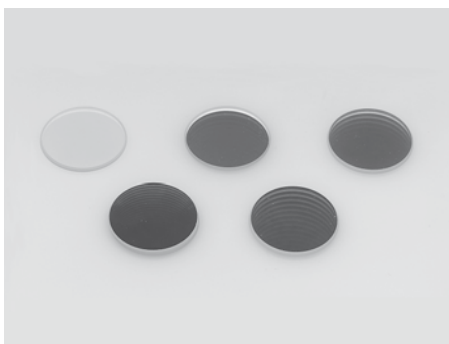

Insert filters for KL 1500 LED and KL 2500 LCD

Insert Filters

- The cold light sources KL 200, KL 1500 LCD, KL 1500 compact, KL 1500 LED and KL 2500 LCD have an integrated filter holder, filter slide or filter wheel.
- Insert filters are placed inside of the light source and evenly illuminate the entire fiber bundle and can be used with all standard fiberoptics for the respective light source.
- Filter change is extremely easy with the filter drawer and filter wheel.
- All SCHOTT insert filters are color interference filters with hard oxide layers, which withstand the extreme temperatures at the focal point of the halogen reflector bulb.
- Manufactured in the EU.

Insert filters				
	Part. No.	Color	Wavelength	Diameter
Insert filter holder KL 200	128 300	As standard integrated in each KL 200		
Insert filters for KL 200	128 302	Blue		Ø 17 mm
	128 303	Red		Ø 17 mm
	128 304	Green		Ø 17 mm
	128 305	Yellow		Ø 17 mm
	128 306	Daylight		Ø 17 mm
Insert filters with holder for KL 1500 LCD/compact	158 302	Blue		
	158 303	Red		
	158 304	Green		
	158 305	Yellow		
	158 306	Daylight		
Insert filter for KL 1500 LED	258 307	Halogene		Ø 28 mm
Insert filters for KL 2500 LCD and KL 1500 LED	258 302	Blue		Ø 28 mm
	258 303	Red		Ø 28 mm
	258 304	Green		Ø 28 mm
	258 305	Yellow		Ø 28 mm
Insert filter for KL 2500 LCD	258 306	Daylight		Ø 28 mm
Fluorescence excitation filters for KL 2500 LCD	258 313	Blue	485 nm	Ø 28 mm
	258 314	Green	515 nm	Ø 28 mm

*Product improvements may result in
specification or feature changes without notice.*

Insert Filters

Lighting and Imaging
SCHOTT AG
Otto-Schott-Strasse 2
55127 Mainz
Germany

Phone: +49 (0)6131/66-7796
Fax: +49 (0)6131/66-7850
lightingimaging@schott.com
www.schott.com/lightingimaging

For your local sales contact please visit
www.schott.com/lightingimaging/english/contact

SCHOTT
glass made of ideas